

Department of Science and Humanities Applied Chemistry Laboratory

Subject: Engineering Chemistry

Observation

Weight of empty crucible=	gm (W1)
Weight of crucible + Sample (Before heating)	= gm (W2)
Weight of sample before drying =	$gm(W_2-W_1)$
=g	gm (W3)
Weight of crucible + sample (after heating) =	gm (W4)
Weight of the sample (after heating) =	gm (W4- W1)
=	gm (W5)
Loss in weight of sample = gm (W5- W3)
= gm ((W6)



Somaiya Vidyavihar University. K. J. Somaiya College of Engineering, Vidyavihar, Mumbai 400077.

Department of Science and Humanities Applied Chemistry Laboratory

Subject: Engineering Chemistry

	Batch:	Roll No.:
	Name:	

<u>Title</u> : Determination of moisture content in Fuel

<u>Aim</u>: To determine the moisture content in the given coal

sample.

Requirement: Porcelain crucible, finely ground charcoal powder,

desiccator, Balance, fractional weight box.

Theory : Moisture is an undesirable component of mined coal. It is

bought and transported at the cost of fuel. It does not contribute to calorific value but actually reduces it. Moisture can be surface moisture which is lost on just drying. However inherent


moisture is not lost by air drying.

Procedure : To determine inherent moisture, air dried coal is crushed

(which can pass through mesh No.60 (ASTM). Initially weigh empty crucible along with lid. Note down the weight. Then weigh about 1 gm of sample in a crucible. Put the lid on the crucible and weigh it again. Note down the weight again. Keep this crucible in an oven maintained at a temperature between 105-110°C. The lid is removed to allow for the evaporation of moisture. Keep it for 1hour in an oven. After 1hour, the crucible is covered with the lid and transferred to a desiccator for cooling. After cooling the crucible is weighed again. Note

down the weight.



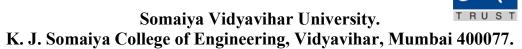

Department of Science and Humanities Applied Chemistry Laboratory

Subject: Engineering Chemistry

<u>Calculat</u>	<u>ions</u>				
Weight of	sample =	gm	(W3)		
Loss in we	eight =		gm (W6)	
		%Moisture	=	Loss in weight Wt of sample take	
			=	<u>W6</u> x 100 W3	
			=		_%
Result	: Percenta	age of moisture	in giv	en charcoal powder	
	_				